


How to implement GS1-traceability in a hospital?





Bernhoven hospital, the Netherlands

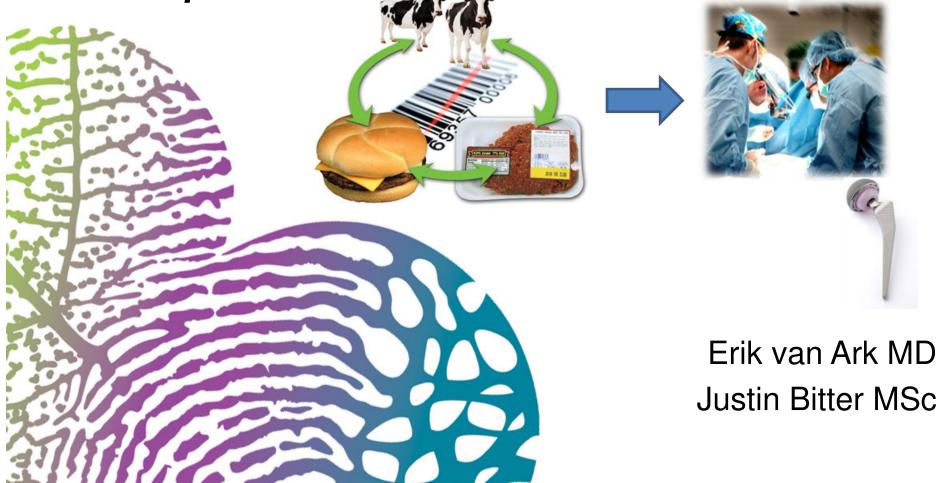




Erik van Ark, MD



- > Anesthesiologist at Bernhoven Hospital, the Netherlands
- Executive director of the medical staff
- Chairman of the OR in a dual management system in which the doctors are in the lead.
- > Innovative: "To Boldly Go Where No Man Has Gone Before"

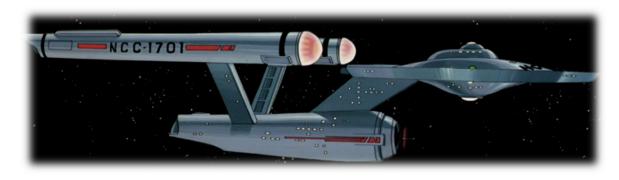

Justin Bitter, MSc




- > Manager OR at Bernhoven Hospital, the Netherlands
- Chairman of the Dutch GS1 focusgroup traceability
- Member of the GS1 Healthcare Leadership Team since July 2014
- > Innovative: "Barcodes are not rocket science"





Baseline Bernhoven..

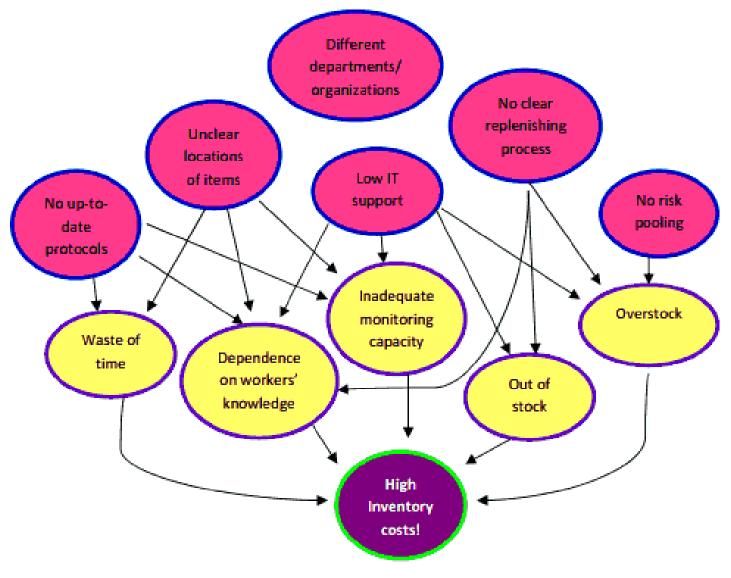




Barcodes







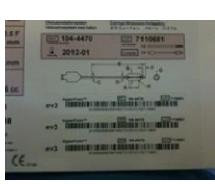




















Traditional thinking

Thinking outside the box



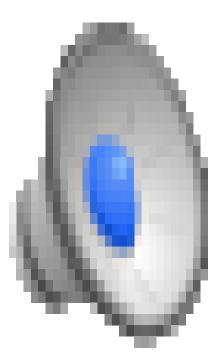



Bernhoven is member of GS1 Netherlands



OR 1 OR 2

OR 3

OR 4

OR 5

OR 6


OR 7

OR 8

OR 9

8718836930195



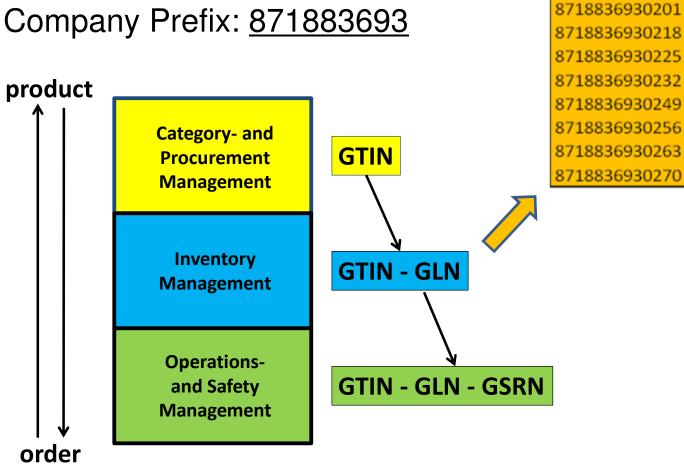



GSTH audit GS1 GO: 10th October 2014




AUDIT INFORMATION					
Objective	Audit the Bernhoven O.R. department to assess compliance to the GS1 Global Traceability Standard for Healthcare (GTSH) utilizing the Global Traceability Conformance (GTC) framework (Issue 1, March 2009)				
Scope	Medical Device Products, specifically implants and sterile disposables. From receipt into OR > OR Logistics process > Patient record; supported by OR Inventory Management System				
Control Points Audited	1.3-1.5, 2.1, 2.3, 3.1-3.4, 3.7, 4.1, 4.3, 4.4, 4.6, 4.7, 4.9-4.11, 5.2-5.5, 5.11-5.14, 6.10-6.11, 6.14, 6.17-6.20, 6.22-6.23, 6.25, 6.27, 6.29-6.30, 7.1-7.6, 8.1-8.4, 9.1-9.3, 10.1-10.3, 10.8, 11.1-11.2, 12.1, 12.3.				
Audit Place	Bernhoven Hospital, Uden, Netherlands				
Audit Date	Friday 10 th October 2014				
Reference documents (Audit Criteria)	GTSH Issue 1.2.0, Oct-2013 Global Traceability Conformance (GTC) framework (Issue 1, March				



Dealing with variation



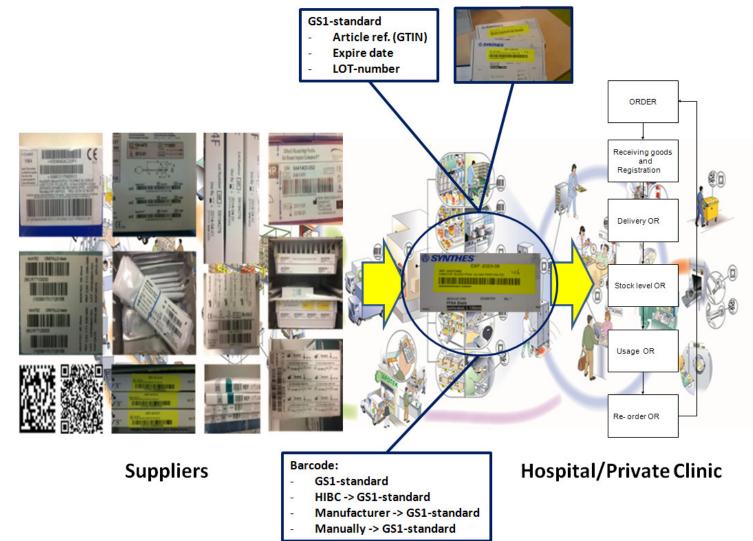



Implants with GS1-128 barcode





Scanning barcodes without relabeling

Bioplug LOT 14-003 Lot de STER. U123	TAILLE SIZE 14	REF H11/A14 2014-06 2019-05	DE PFNA - Pr ES PFNA - Cl FR PFNA - Cl IT PFNA - Ch	S 9140816 ox. Fem. Nail Ø 10mm, L 240mm ox. Femurnagel Ø 10mm, L 240mm avo prox. fem. Ø 10mm, L 240mm ou prox. p/fémur Ø 10mm, L 240mm hiodo fem. pross. Ø 10mm, L 240mm wilha prox. fem. Ø 10mm, C 240mm	Qty. 1 Date of sterilization: 2014-09 2024-09 STERILE R
(01)03770000362	2039(17)1905	01 (10)14-003U123	Manufacturer Synthes GmbH Eimattstrasse 3 4436 Oberforf Swtzerland	(01) 07611819289239 (17) 240901 (10) 9140816	Material: Ti Al6 Nb7



Implants with damaged barcode



Relabeling barcode with GTIN of company



Implants not GS1-128 compliant


REF: 12763026, HEUPSTEEL SP CCD 126 L3 90701 LOTNO: 1405203523 140520/3 Relabeling to GS1-128 SN barcode by generating own REE 127-630/26 **GTIN** 2763026, HEUPSTEEL SP CCD 126 La 190701 LOTNO: 1405203523

11:58



Call for standardization of barcodes



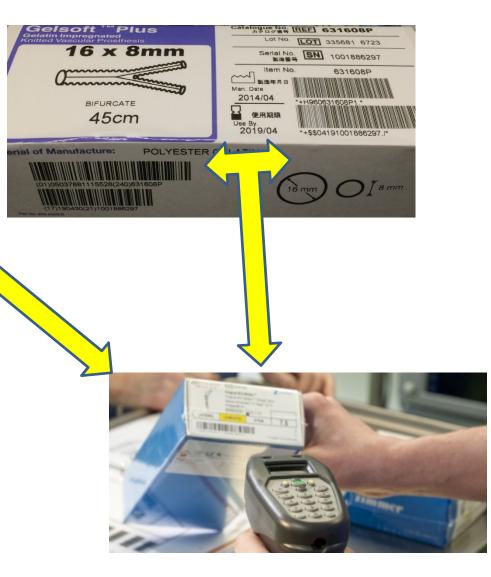





Bar Code Symbology

The familiar symbology (format) of bar codes, particularly in the retail environment, is that of the '1D' or 'linear' type, vertical black lines and white spaces (Figure 5), but over recent years another generation of bar code symbology is being increasingly used, the "two-dimensional (2D) or "matrix" type symbology, e.g. GS1 DataMatrix (Figure 3).



Figure 5




These two symbologies contain the same data.

Figure 6

http://www.gs1.org/sites/default/files/docs/healthcare/20130725%2 0HPAC%20Position%20Paper%20Bar%20Code%20Issues.pdf



Nationwide register for implants

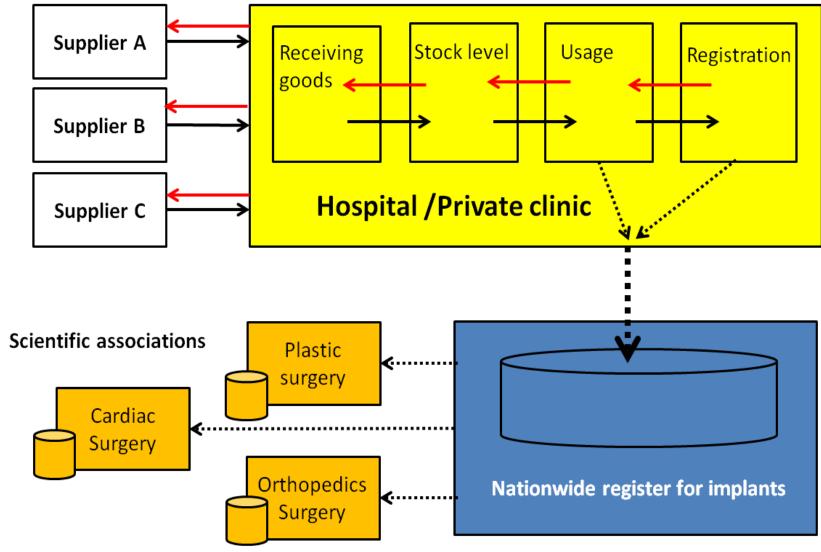



Nationwide register for implants

- (8018)871883693064626937(8008)140801150100(01)05050474558199(17)190300(10)M1203005
- (8018)871883693064626937 Company Prefix/Patient
- (8008)140801150100 Year/Month/Day/Hour/Minute/Second
- (01)05050474558199 Global Trade Item Number
 - (17)190300 Expire date
- (10)M1203005 Lot number



GS1 focusgroup traceability in Healthcare



Agreements on the use of GS1 standards



Next steps (2015) in Bernhoven Hospital:



High risk medication

Instrument trays

Thank you for your attention!!

21

